
uSqueez 3

OS-3208 Operation Manual

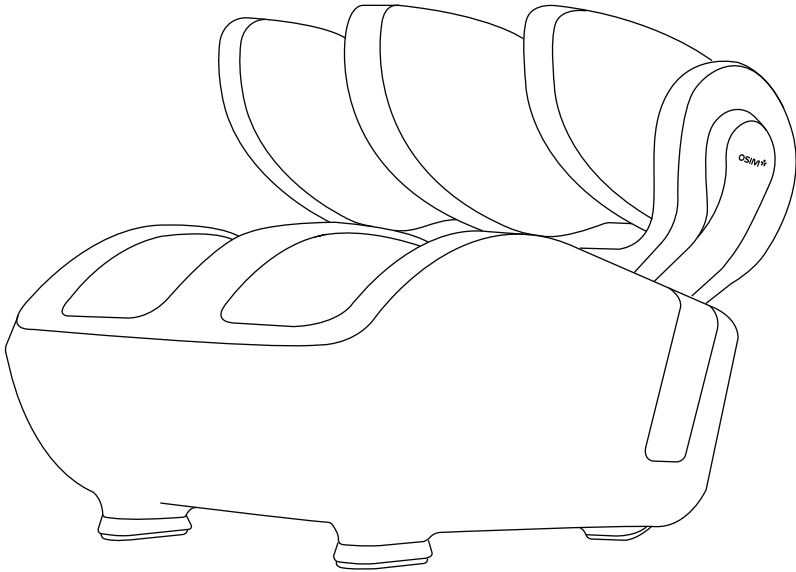1) English

2) Simplified Chinese 简体中文

3) Traditional Chinese 繁體中文

OS-3208 uSqueez 3

Thank you for purchasing the OSIM uSqueez 3. Please read the operation manual thoroughly for the correct methods of usage and take special notice of the safety section before you start.

This manual should be kept on hand for easy reference. It is your guide to a safe and efficient operation.

Contents	Safety	2
	Names of parts and functions	5
	Preparation before use	7
	Operation	8
	Product care and storage	14
	Troubleshooting	15
	Specifications	16

IMPORTANT SAFETY INSTRUCTIONS

When using the product, please **take special notice** of the following basic precautions:
Read all instructions before using this product.

DANGER – To reduce the risk of electric shock:

Always unplug this product from the electrical outlet immediately after use and before cleaning.

WARNING – To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1) Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- 2) The product should never be left unattended when plugged in. Unplug from electrical outlet/power source when it is not in use, and before putting on or taking off parts.
- 3) Do not operate the product under the blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 4) Do not use the product if the electrical outlet/power source is loose.
- 5) Do not connect or disconnect the power plug/USB cable from the electrical outlet/power source with wet hands.
- 6) When disconnecting the power plug/USB cable, pull it out by holding the plug/USB cable head, not the cord/cable.
- 7) Do not operate the product if it has a damaged cord or plug /USB cable, a malfunction, a power failure, or it is dropped into the water. Injury may occur when power is restored unexpectedly. Return the product to OSIM Service Center for examination and repair.
- 8) If the power plug/USB cable or power cord is damaged, it must be replaced by OSIM Service Center or its authorized service agents.
- 9) Do not carry this product by the power cord /USB cable or use the power cord/USB cable as a handle.
- 10) Keep the product power cord/USB cable away from heated surfaces.
- 11) Do not place the power cord/USB cable under the product or any other heavy objects.
- 12) Do not wind the power cord/USB cable around the product.
- 13) Do not operate the product with the air openings blocked. Keep the air openings free of lint, hair, and the likes.
- 14) Remove power plug/USB cable from the electrical outlet /power source immediately should liquid be accidentally spilled onto the product as this may result in electric shock.
- 15) Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- 16) To disconnect, turn off all controls and remove the power plug/USB cable from the electrical outlet/power source.
- 17) Do not drop or insert any object or flammable material into any openings of the product. This may result in electric shock, malfunction or risk of fire.

18) Children should be supervised to ensure that they do not play with the product.

19) This product can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Keep children away when the product is in operation.

20) Cleaning and user maintenance of the product shall not be made by children.

21) Do not use or operate the product when your body or hands are wet.

22) Do not use this product outdoors.

23) Do not use the product in an excessively humid or dusty environment as it may result in malfunction.

24) Do not use the product in temperatures higher than 40°C.

25) Do not expose the product to heaters, stoves or direct sunlight.

26) Do not cover the product during operation. This may result in a rise in motor temperature and subsequent malfunction or fire.

27) Do not use the product in the event of a lightning storm.

28) If you should use the product in a cold room, do not increase the room temperature abruptly. It is recommended that you increase the temperature gradually to a normal level.

29) If the product is brought into a warm environment from a cold place where it was stored, wait for an hour before use as there may be water droplets on mechanical parts due to condensation. Using the product under such condition may result in malfunction.

30) When first using the product, please be careful and test it out using a low intensity/or use within a short period of time to understand how the different motors work and your preferred intensity.

31) Be sure to start on the gentle massage program if using the product for the first time.

32) If your legs slip off the footrest during massage, do not force them back in when the kneading discs are clamped as this may result in injury.

33) Do not stand on or in product. Use only while seated.

34) Place and use the product on an even, non-slip floor.

35) Do not use the product if the fabric of the footrest is torn and the internal components are exposed.

36) Do not disassemble the fabric cover of the product. Touching the internal components of the product may result in malfunction or an electric shock.

37) Do not dismantle, drop, or subject the product to shock or vibration.

38) Do not use the product in combination with other therapeutic equipment or electric blanket, etc, as this may result in ineffectiveness or injury.

39) Do not rest or place heavy objects on the product as this may result in malfunction.

40) The product must only be serviced or repaired by the OSIM Service Center. Users must not attempt to disassemble or repair the product.

41) Some noise may be heard coming from the product while in use. This is due to the structure of the product and is normal.

- 42) Please do not use the product if you are not feeling well or under the influence of alcohol.
- 43) Stop using immediately should you experience discomfort and consult your doctor if you experience any abnormal or severe pain.
- 44) This product is not intended for self-diagnosis or self-treatment of conditions that should be managed by a qualified health care provider.**
- 45) If you are under any medication or fall into any of the following categories (but not limited to), please consult your health care provider first before purchase use:**
- Heart diseases
 - High, or low blood pressure, deep vein thrombosis (DVT)
 - Pregnant
 - Diabetes
 - Individuals fitted with medical implants (e.g. pacemaker, automatic cardiac defibrillator (AICD))
 - Recently underwent surgical operations
 - Those with any bone conditions like fractures, osteoporosis
- 46) The suggested usage duration is 15 minutes. Ensure that the product rests for 10 minutes before it is being used again as continuous usage over a long period may activate the overheat protection device.
- 47) It is recommended that your usage do not exceed a total of 30 minutes per day as this may result in excessive stimulation and adverse effects.
- 48) The product has a heated surface. People who are insensitive to heat must be careful when using the product.



SAVE THESE INSTRUCTIONS

This product is intended for household use only

Warranty is subjected to cancellation if the product is used for commercial purposes.

Names of parts and functions

Main unit

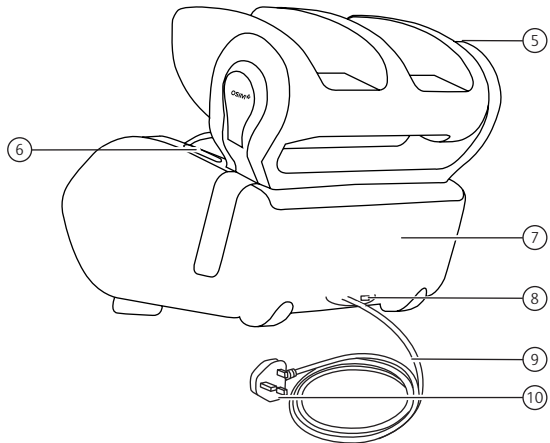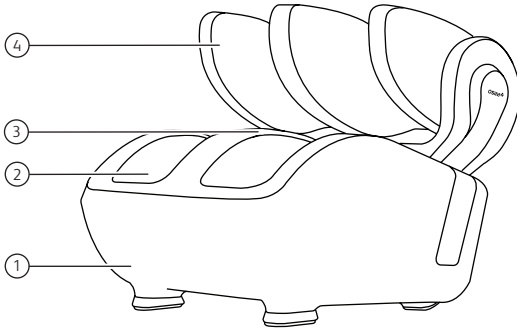*Product images are for illustrative purposes only and are not drawn to scale.


1. Bottom base
2. Washable fabric lining with zip
3. Front casing
4. Kneading plate
5. Control panel
6. Safety switch (safety switch will be activated once there is object or external force applied on it, then kneading plate will stop moving downwards and rotating down)

7. Back casing
8. Main switch
9. Power cord
10. Power plug*

*Power plug will vary according to country

Names of parts and functions

Control panel

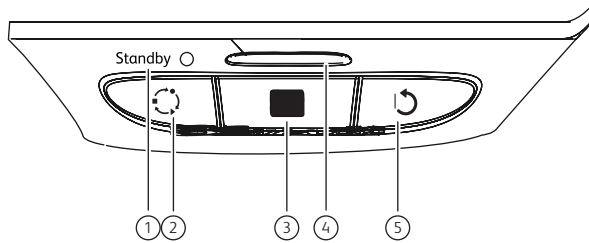
1. Standby LED indicator
2. [🔄] 'Auto Program' button (Three auto programs 'Japanese Shiatsu', 'Thai Stretch' and 'Chinese Tui-Na' can be selected through control panel. To enjoy more auto programs, please control via App.)
3. [■] 'Stop' button
4. Power LED indicator
5. [🔄] 'Reset' button

Preparation before use

Moving the uSqueeze 3

- Lift the uSqueeze 3 by holding onto the two sides of the product as shown.
- The uSqueeze 3 must be completely lifted off the floor when moving from one place to another. Do not drag the uSqueeze 3.
- Practice good lifting body posture when moving the uSqueeze 3 to protect your back.


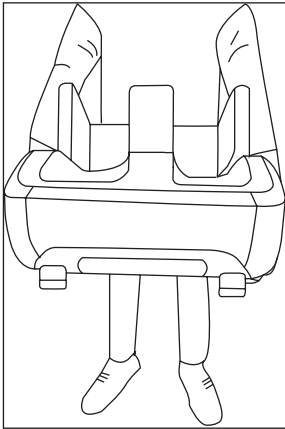
Caution

- Ensure that the power cord is not placed under the uSqueeze 3.
- Always place the uSqueeze 3 on even, non-slip floor. For floors that may be easily scratched, you may place a mat under the uSqueeze 3.

Operation

To switch on the power

- Connect the power plug to the electrical outlet.
- Turn on the main switch located at the lower back of the uSqueeze 3.
- The Power LED indicator on the control panel will be illuminated.

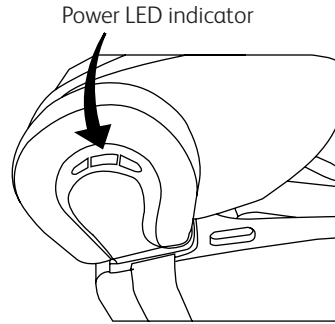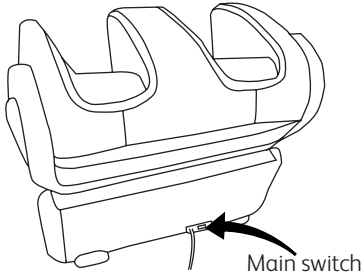
Control Panel

Explaining the Power LED light indicator

Power	Red Light
Bluetooth Pairing	Blue Light
Auto Program Japanese Shiatsu	White Light
Auto Program Thai Stretch	Yellow Light
Auto Program Chinese Tui-Na	Green light

To start the auto massage program on control panel

- Place your legs into the uSqueeze 3.
- Press the [🔄] button once to start 'Japanese Shiatsu' program (at Wide calf setting).
- Press the [🔄] button again to start 'Thai Stretch' program (at Wide calf setting).
- Press the [🔄] button again to start 'Chinese Tui-Na' program (at Wide calf setting).
- Press the 'Stop' button [■] to activate emergency stop for the auto massage programs at anytime during operation.
- Press the 'Reset' button [↺] to normally stop the auto massage programs at anytime during operation and kneading plate will go home position.
- Please download OSIM Relax & Relieve App to explore more functionality on uSqueeze 3.

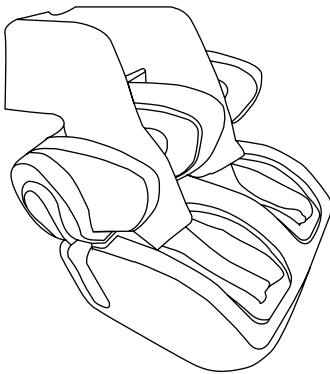Operation

Pairing iOS and Android devices with OSIM uSqueeze 3

Relax & Relieve

Ensure Bluetooth is enabled on the uSqueeze 3



- Ensure that the main switch on the OSIM uSqueeze 3 are turned on.
- Bluetooth will be enabled once the main switch is turned on.

Pairing iOS and Android devices with OSIM uSqueeze 3

Initiate Bluetooth pairing on your iOS and Android device.

- Turn on the Bluetooth function on your smart device.
- Tap on the OSIM Relax & Relieve App on your smart device.
- Tap on 'Leg massager'. The Relax & Relieve App will search for the unique device ID.
- Select the correct Device ID as shown, select 'OS3208 XXXXXXXX' (XXXXXXX can be any numerical values) to connect with the uSqueeze 3.

Operation

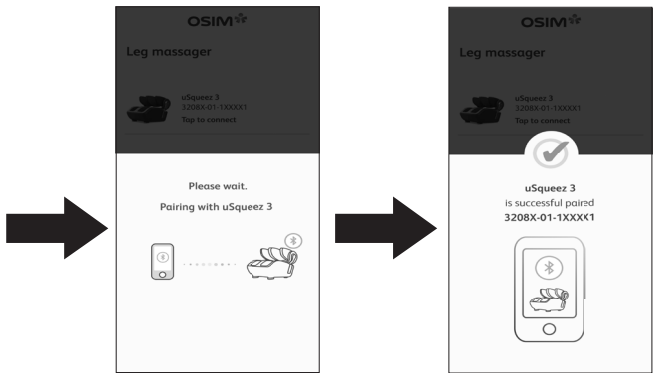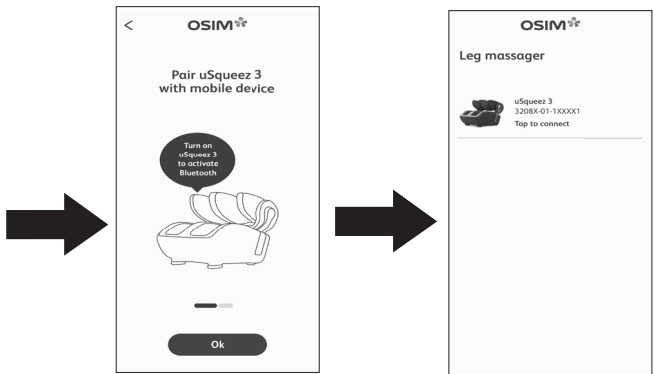Operation

Download new massage programs to OSIM uSqueeze 3

1. After you have successfully paired your smart device with uSqueeze 3, tap on the Downloads icon [📁] in the main menu.
2. At the Downloads screen, a list of new massage programs will be displayed.

Notes

- iOS or Android device needs to have network connection for new massage programs to be downloaded.
- 3. Tap on the desired massage program to preview the massage program details.
- 4. Tap on 'Download' to download the free massage program or 'Cancel' to return to previous screen.

Functionality of OSIM Relax & Relieve App

User is able to access the following through the OSIM Relax & Relieve App.

- New downloadable massage programs will be pushed to App as Downloads when available to keep your OSIM uSqueeze 3 always updated with latest massage programs.

Compatibility of OSIM Relax & Relieve App

- Compatible with selected smart devices running on iOS version 11 and above or Android version 8.0 and above.
- Due to variations in operating systems, non-compatibility could occur. We recommend updating your phone's operating system whenever a new version is available to enjoy the new App features. The minimum phone operating system requirement are subject to change without notice.
- All smart devices require a compatible wireless connectivity (Wi-Fi or cellular data networks). During setup process, make sure to connect your mobile device to 2.4 GHz network as the devices are only compatible with 2.4 GHz networks. Once setup completed, you can continue to operate device using your regular Wi-Fi or cellular data networks.
- For updates on OSIM Relax & Relieve App, please visit the respective app stores.
- Notifications of app updates should not be disabled, and it is recommended that the user installs all updates to ensure that the app will function correctly and is kept up to date.

Operation

- In the event of power failure, switch off the main switch and disconnect the power cord to prevent any sudden operation if the power is restored unexpectedly.
- Do not stand or exert excessive pressure on the uSqueeze 3. This may significantly slow down or stop the motor, potentially causing overheating and subsequent damage to the uSqueeze 3.
- There is an Auto-Timer of 15 minutes for all massage modes.
- After each 15 minutes usage, ensure the uSqueeze 3 rest for at least 10 minutes before starting the next session. If uSqueeze 3 is operated consecutively for a period of time, it may activate the overheat protection device and stop operating automatically.
- If uSqueeze 3 stops operating due to over heating, wait for the temperature to cool down approximately 30 minutes before you begin again.
- Press the Stop button [■] to terminate the massage if you experience any discomfort. The massage will stop immediately, kneading discs will complete its last cycle before returning to the storage position.
- Product comes with a warm air function and can be activated via the OSIM Relax & Relieve App.

Product care and storage

Product placement and storage

Environmental factors can strongly affect the product physically and mechanically. Do not expose the product to direct sunlight or in places with extreme heat or air-conditioning vents. These can cause the product color to fade.

If the product is not to be used for a long period of time, cover it with a dust cover and disconnect the power plug from the electrical outlet.

Clean the product using a dry cloth regularly

Before cleaning, ensure that the main switch is turned off and the power plug is disconnected from the electrical outlet.

The product is customarily exposed to the natural body oils of users. The main unit structure and fabric can be cleaned thoroughly with a piece of soft, dry cloth. The common microfiber cleaning materials or gloves are very well suited with this routine maintenance process. Such materials and gloves are available for sale in most DIY shops.

Do NOT use wet elements for cleaning

The product is made up of electrical components and mechanical parts, which are all very sensitive to moisture. Under no circumstances should they be exposed to any form of water or moisture. If you use water (eg. a wet cloth) to clean the product, the coating may be eroded and the protection for the product will be diminished.

As for the control panel, it may result in malfunction if you use a wet cloth to clean it.

Do NOT use cleaning solvents

Under no circumstances should you use any cleaning agents, like benzene or thinner to clean the product. This may cause the color of product to fade.

Do not wash the fabric using domestic washing machine. This may result in the deterioration of quality.

Do NOT attempt to repair your product

Under no circumstances should you attempt to repair the product, for safety reason. If you suspect a malfunction, always call the OSIM Service Center for further assistance and information.

Troubleshooting

If you encounter problems with the operation of the uSqueeze 3, please check the following points to see if they can be rectified. If the uSqueeze 3 still malfunctions or does not work, please turn off the uSqueeze 3 and remove the power plug from the electrical outlet. Contact OSIM Service Center* for servicing.

Do not dismantle or attempt to repair the uSqueeze 3 yourself. Doing so will void the warranty and OSIM will not bear any responsibility for any subsequent consequences.

Problem	Possible causes / action
<ul style="list-style-type: none">• Motor sound during operation.	<ul style="list-style-type: none">• The sound is the result of the mechanical structure in operation and not a malfunction.
<ul style="list-style-type: none">• The massage stops midway.	<ul style="list-style-type: none">• The Auto Timer function of 15 minutes may have lapsed.• The 'Stop' button [■] on the control panel or in the App was pressed accidentally.• The overheat protection device may have been activated. Turn off the main switch and allow the uSqueeze 3 to cool down at least 30 minutes before using it again.
<ul style="list-style-type: none">• The uSqueeze 3 will not start.	<ul style="list-style-type: none">• Ensure that the power plug is firmly connected to the electrical outlet.• Ensure that the main switch is turned on.

* Please visit us at www.OSIM.com for information on the nearest service center. Product warranty is valid only in the country of purchase.

Specifications

Model name	OSIM uSqueez 3
Model number	OS-3208
Voltage available	220-240V 50-60Hz (95W) 110V ~ 60Hz (95W) 220V ~ 50Hz (310W) <small>(Kindly ensure that the local Alternating Current (AC) voltage corresponds with product specified voltage)</small>
Auto timer	15 minutes
Dimension	54 cm (L) x 54.5 cm (W) x 48 cm (H)
Net weight	Approx. 16 kg
Gross weight	Approx. 19 kg

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