
uErgoh

OS-8245 Operation Manual

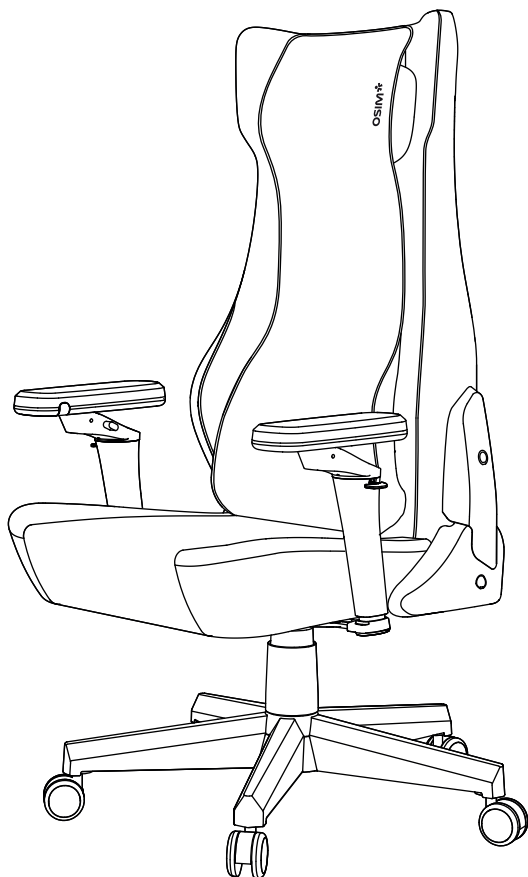1) English

2) Simplified Chinese 简体中文

3) Traditional Chinese 繁體中文

4) Specification 规格说明

OS-8245 uErgoh

Thank you for purchasing the OSIM uErgoh. Please read the operation manual thoroughly for the correct methods of usage and take special notice of the safety section before you start.

This manual should be kept on hand for easy reference. It is your guide to a safe and efficient operation.

Contents	Safety	2
	Names of parts and functions	5
	Preparation before use	7
	Operation	13
	App download	15
	Product care and storage	17
	Troubleshooting	18

IMPORTANT SAFETY INSTRUCTIONS

When using the product, please take special notice of the following basic precautions:
Read all instructions before using this product.

DANGER – To reduce the risk of electric shock:

Always unplug this product from the electrical outlet immediately after use and before cleaning.

WARNING – To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1) Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- 2) The product should never be left unattended when plugged in. Unplug from electrical outlet/power source when it is not in use, and before putting on or taking off parts.
- 3) Do not operate the product under the blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 4) Do not use the product if the electrical outlet/power source is loose.
- 5) Do not connect or disconnect the power plug/USB cable from the electrical outlet/power source with wet hands.
- 6) When disconnecting the power plug/USB cable, pull it out by holding the plug/USB cable head, not the cord/cable.
- 7) Do not operate the product if it has a damaged cord or plug/USB cable, a malfunction, a power failure, or it is dropped into the water. Injury may occur when power is restored unexpectedly. Return the product to OSIM Service Center for examination and repair.
- 8) If the power plug/USB cable or power cord is damaged, it must be replaced by OSIM Service Center or its authorized service agents.
- 9) Do not carry this product by the power cord/USB cable or use the power cord/USB cable as a handle.
- 10) Keep the product power cord/USB cable away from heated surfaces.
- 11) Do not place the power cord/USB cable under the product or any other heavy objects.
- 12) Do not wind the power cord/USB cable around the product.
- 13) Do not operate the product with the air openings blocked. Keep the air openings free of lint, hair, and the likes.
- 14) Remove power plug/USB cable from the electrical outlet /power source immediately should liquid be accidentally spilled onto the product as this may result in electric shock.
- 15) Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- 16) To disconnect, turn off all controls and remove the power plug/USB cable from the electrical outlet/power source.
- 17) Do not drop or insert any object or flammable material into any openings of the product. This may result in electric shock, malfunction or risk of fire.
- 18) Children should be supervised to ensure that they do not play with the product.
- 19) This product can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Keep children away when the product is in operation.
- 20) Cleaning and user maintenance of the product shall not be made by children.

- 21) Do not use or operate the product when your body or hands are wet.
- 22) Do not use this product outdoors.
- 23) Do not use the product in an excessively humid or dusty environment as it may result in malfunction.
- 24) Do not use the product in temperatures higher than 40°C.
- 25) Do not expose the product to heaters, stoves or direct sunlight.
- 26) Do not cover the product during operation. This may result in a rise in motor temperature and subsequent malfunction or fire.
- 27) Do not use the product in the event of a lightning storm.
- 28) If you should use the product in a cold room, do not increase the room temperature abruptly. It is recommended that you increase the temperature gradually to a normal level.
- 29) If the product is brought into a warm environment from a cold place where it was stored, wait for an hour before use as there may be water droplets on mechanical parts due to condensation. Using the product under such condition may result in malfunction.
- 30) When first using the product, please be careful and test it out using a low intensity/or use within a short period of time to understand how the different motors work and your preferred intensity.
- 31) Be sure to start on the gentle massage program if using the product for the first time.
- 32) Do not stand on or in product. Use only while seated.
- 33) Place and use the product on an even, non-slip floor.
- 34) Do not use the product if the fabric of the footrest is torn and the internal components are exposed.
- 35) Do not disassemble the fabric cover of the product. Touching the internal components of the product may result in malfunction or an electric shock.
- 36) Do not dismantle, drop, or subject the product to shock or vibration.
- 37) Do not use the product in combination with other therapeutic equipment or electric blanket, etc, as this may result in ineffectiveness or injury.
- 38) Do not rest or place heavy objects on the product as this may result in malfunction.
- 39) The product must only be serviced or repaired by the OSIM Service Center. Users must not attempt to disassemble or repair the product.
- 40) Some noise may be heard coming from the product while in use. This is due to the structure of the product and is normal.
- 41) Please do not use the product if you are not feeling well or under the influence of alcohol.
- 42) Avoid food intake for at least an hour before massage /usage.
- 43) Each auto massage program is designed to last 15 -30 minutes for optimal effectiveness. Ensure that the product rests for 10 minutes before it is being used again as continuous usage over a long period may activate the overheat protection device.
- 44) For best results, enjoy your massage regularly, while keeping in mind that moderation is key to avoid excessive stimulation and ensure a balanced experience.
- 45) Stop using immediately should you experience discomfort and consult your doctor if you experience any abnormal or severe pain.
- 46) Do not exert excessive force on the back of your neck while using massager.
- 47) This product is not intended for self-diagnosis or self-treatment of conditions that should be managed by a qualified health care provider.

- 48) If you are under any medication or fall into any of the following categories (but not limited to), please consult your health care provider first before use:
- Heart diseases
 - High, or low blood pressure, deep vein thrombosis (DVT)
 - Pregnant
 - Diabetes
 - Individuals fitted with medical implants e.g. pacemaker, automatic implantable cardioverter defibrillator (AICD)
 - Recently underwent surgical operations
 - Those with any bone conditions like fractures, osteoporosis
- 49) This product contains lithium battery that is only replaceable by a skilled technician. Do not attempt to assemble or repair the product. Contact OSIM Service Centre for any repair.
- 50) The battery pack must be removed from the product before discarding - the product must be powered off before removing the battery pack. Return the battery pack to OSIM Service Center for recycling or disposed of in accordance with local waste management regulations.
- 51) Do not overcharge otherwise it may damage the battery and affect usage.
- 52) Never use power strips, or extension cords to charge your lithium batteries.
- 53) Stop charging once the product is fully charged.
- 54) Keep product away from items that can catch fire, including bedding and couches. Never charge on surfaces that can catch fire.
- 55) Never leave it plugged in when it is not charging.
- 56) Never charge the product near doors or hallways that may block exits.
- 57) Do not charge or store your lithium-ion battery in extreme temperatures.
- 58) Handle your lithium-ion battery with care. Do not drop or submerge in water.
- 59) Do not modify or attempt to repair your lithium-ion battery.
- 60) Stop using the battery if you notice these problems:
- Burning electrical odor emitted from the battery
 - Change in the battery or product shape
 - Change in the battery color
 - Noises or leaking coming from the battery


SAVE THESE INSTRUCTIONS

This product is intended for household use only

Warranty is subjected to cancellation if the product is used for commercial purposes.

Names of parts and functions

Main unit

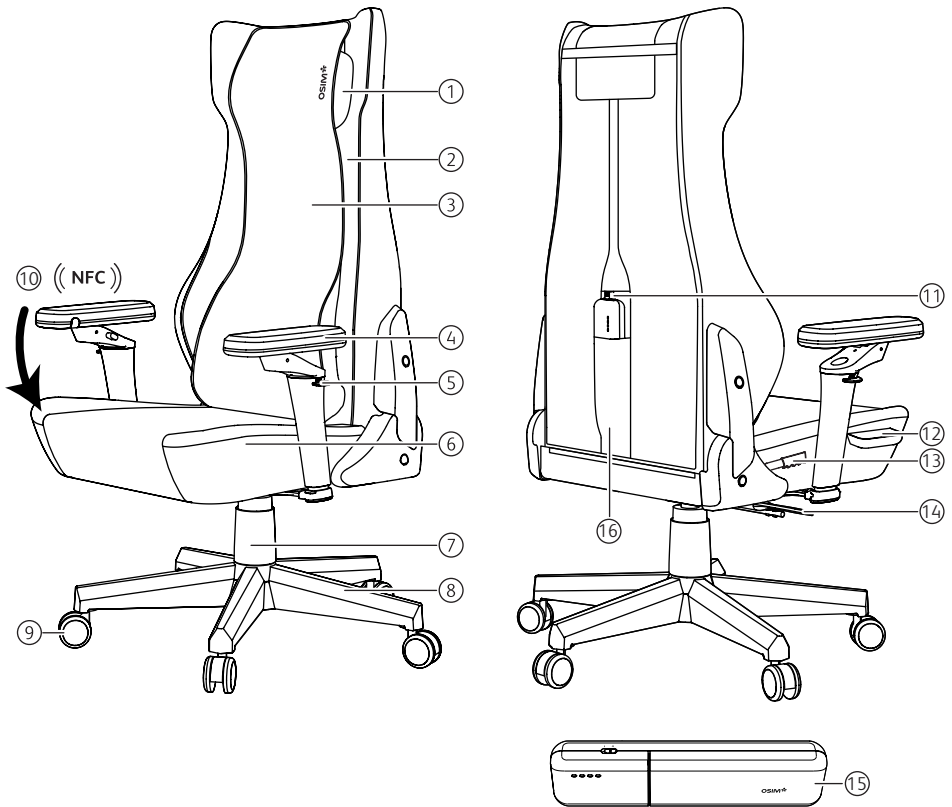
*Product images are for illustrative purposes only and are not drawn to scale.

- | | |
|---------------------------------|-------------------------------|
| 1. Head cushion | 10. NFC |
| 2. Backrest | 11. Power source connector |
| 3. Backrest pad | 12. Control panel |
| 4. Armrest | 13. Backrest locking lever |
| 5. Armrest height locking latch | 14. Seat height locking lever |
| 6. Seat | 15. Battery pack* |
| 7. Hydraulic piston | 16. Battery holder |
| 8. Wheelbase | |
| 9. Caster wheels | |

* Only use the provided USB-C cable for charging the battery pack.

Names of parts and functions

Control Panel

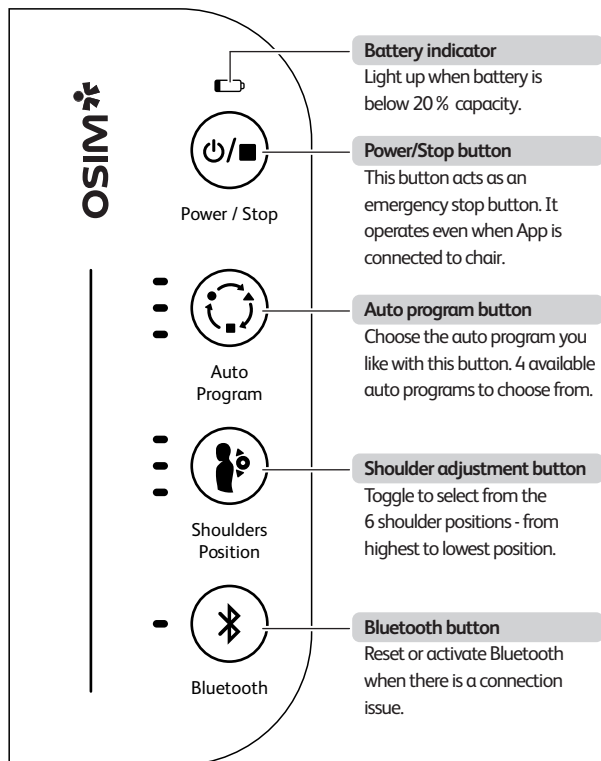
Caution

- Do not sit, step or place heavy objects on the control panel.
- Do not exert excessive force on the control panel.
- Do not operate the control panel with wet hands.
- Do not operate the control panel using sharp-edged or hard objects. It may result in damage.

Preparation before use


Setting up the product

- Unpack the carton and assemble the items accordingly to the provided assembly guide.
- Ensure a space allowance of at least 80cm behind the product to prevent obstruction when it reclines.

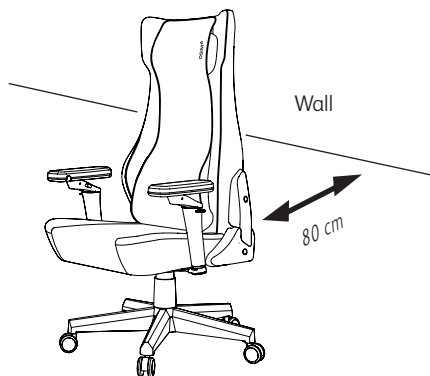Placing the backrest pad

- Placing the backrest pad on the seat as shown, insert its flap through the strap in the back.

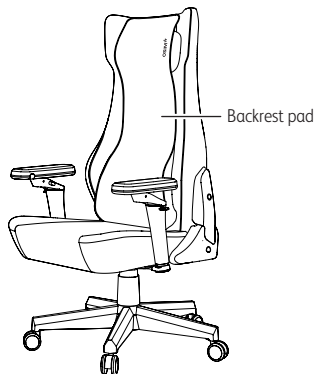
Note

- Remove backrest pad before starting massage.

Placing the head cushion


- Place the head cushion on top of the backrest and insert its flap through the strap as shown.

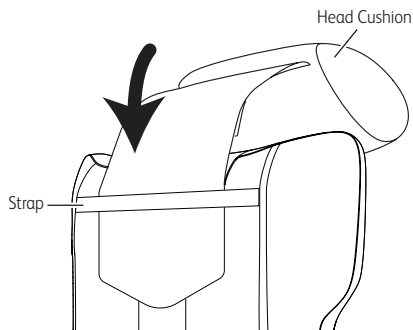Preparation before use

Charging the battery pack

- Before first use, it is recommended to fully charge the battery pack after unboxing, for better user experience.
- When battery is removed from the packaging, the Power On/Off switch will be set to [**O**]. Slide the switch to [**I**] (the battery LED indicator will light up) before charging it with the provided USB-C cable.
- It takes approximately 4 hours to fully charge a new battery. Using a USB Fast Charging adapter, the battery takes about 4 hours to fully charge; otherwise, it may take up to 8 hours.
- When fully charged, the battery will display four white LEDs, which will turn off after approximately one minute.
- Each full charge last at least 60 minutes of usage on average (Usage time varies based on the massage program, intensity and environmental temperature).

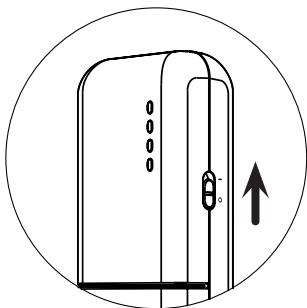
Caution

- The product is inoperable during charging.
- Please remove the battery pack from the holder before charging.
- Please only use the provided USB-C cable for charging.
- Please do not leave the charging battery pack unattended - keep it away from flammable materials.
- Allow the battery pack to cool down before and after charging.
- Please unplug the battery pack from the USB-C cable after it is fully charged.
- Ensure the power cable is plugged to the battery pack in the correct direction - check the orientation of the battery pack socket if it takes much force to plug it in.
- The product is to be operated within the recommended environmental temperature.
- During charging, if excess heat is noticed on the battery pack or USB-C cable, stop charging immediately. Unplug the USB-C cable from battery pack and do not keep the battery pack in the holder. Contact OSIM Service Center immediately for advice.
- The battery pack is only serviceable by OSIM Service Center - do not attempt to open it.

Preparation before use

Keeping the battery pack in the battery holder at the back of product

- Plug the power cable (above the battery holder) to the battery pack.
- Take note of the orientation of the battery pack relative to the power cable. The power cable will only fit into the battery pack when oriented correctly. (The battery LED indicator is located at the top and front facing).

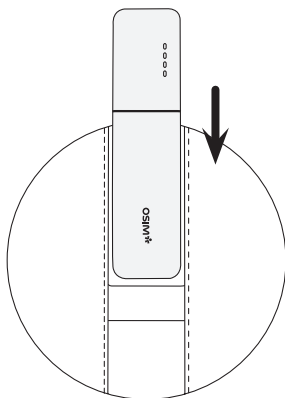
Note

- Do not place heavy objects on the battery.
- If the product will not be in use for 2 months or more, remove the battery from the product. Slide the Power On/Off switch to [**O**] before storage. Please note that battery pack should be charged at least once every 6 months even if it is not in use, to maintain their longevity.
- Before using the product again, slide the Power On/Off switch to [**I**] on the battery and fully charge it before using it with product.

Moving the product

- The product is fitted with casters for easy of manoeuvring over a short distance.
- Raise the backrest fully to its upright position and lower the seat height.
- Push the product carefully (shown below) to the desired location.

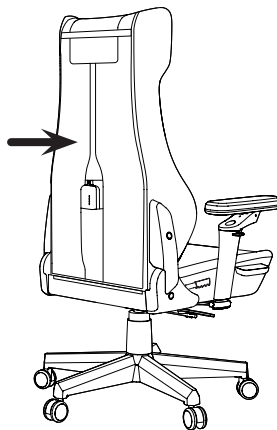
Note

- Do not move the product with someone on it as this might tip over the product and cause an accident and injury.

Preparation before use


Seat height adjustment

- Push down the seat height locking lever to adjust the seat height.
- Without any force of gravity, the seat height will rise to its maximum. If seated in the product, the seat height will be lowered gradually.
- Push down the lever without any weight on the seat to rise the seat height to the maximum. Push down the lever when seated, release the lever to lock upon reaching the desired seat height.

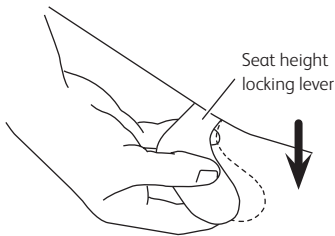Backrest angle adjustment

- To recline the backrest, pull the backrest locking lever upward.
- Lean back slowly to avoid possible toppling.
- Be mindful when raising the backrest, as it may spring back to upright position - Lean forward slowly to prevent the back or head from being hit by the massage mechanism.

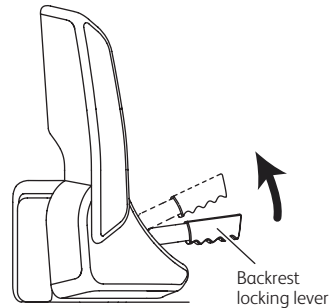
Caution

- Ensure that the backrest is at its most upright position before sitting in the chair and reclining to your preferred position.
- Do not place fingers or objects in the gap between the backrest and seat during adjustment.
- Ensure that there is no obstruction when reclining the backrest.


Preparation before use

Armrest position adjustment



- Lift the latch to adjust the height of the armrest.

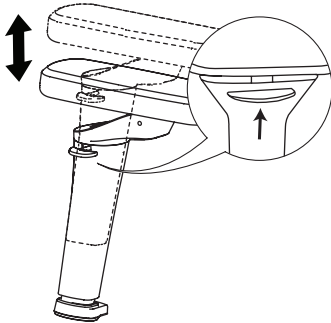
- Press the bottom button to tilt up the armrest angle.

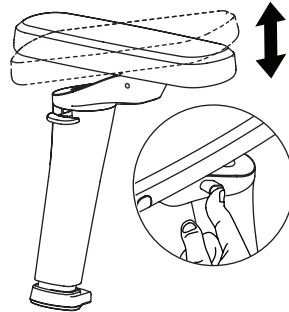- Press the top button to adjust the armrest forward/backward and swivel 360°.

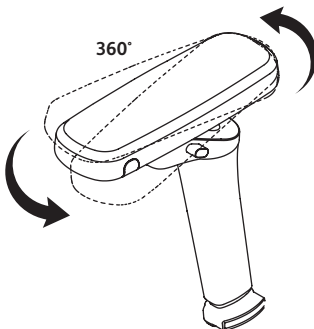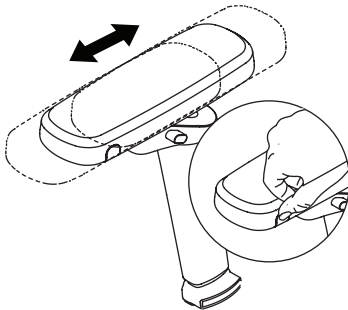
- The armrest base can be rotated freely by hand for a wider positioning.

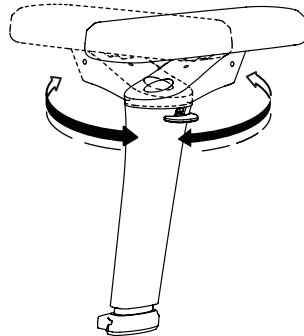
Preparation before use

Usage and safety guidelines

Improper usage may damage the product or cause serious injuries. To ensure your safety while using the product, adopt the following guidelines:

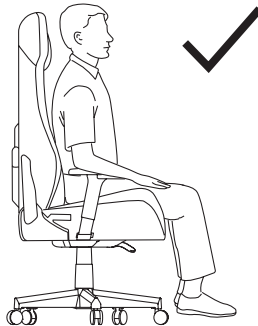- Always sit in the center of the product, with your back against the backrest.

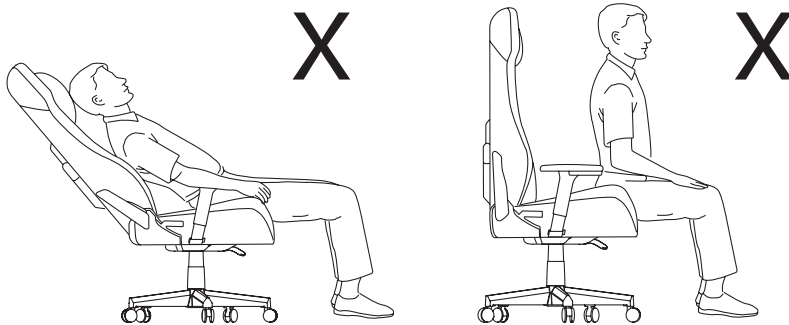- Do not sit on the front edge of the product, this may cause chair to tip over and result in injury.

Caution

- This product is designed for one person at any one time.
- Do not stand or squat on the product.
- Do not sit on the armrests.
- Do not use the product to move other objects.
- Ensure product is used on a flat floor surface.
- Ensure product is properly assembled using the provided assembly guide.
- Check that all parts are in good condition regularly.
- Do not use the chair if there are any missing, damaged, or worn components.

Operation

Chair control via control panel

After plugging in the battery, the control panel lits up and the product is on Standby mode. If there is no operation/App connection detected within 4 minutes, it will enter Sleep mode.

- To exit Sleep mode
 - Press the Power/Stop [ - Together with the battery indicator, both the Shoulders Position and Bluetooth LEDs will lit up.
 - The Shoulder Position LED blinks as the product resets.
- To operate the product
 - Press the “Shoulders Position” [ - Once the correct shoulder position is set, press the “Auto Program” [ - You are all set to relax and enjoy your massage.
- If operating the product with the control panel, the Bluetooth will go into Sleep mode after 4 minutes of inactivity. In any time when user wants to operate the product via OSIM Well-Being App, simply press and hold the Bluetooth button [- The product has a “low battery warning” function.
 - When the battery capacity falls below 20 %, the battery indicator LED will light up - with a long beep.
 - When no action is taken, another long beep will be sounded about a minute after the first beep.
 - If there is no further action, the product may shut down halfway through a massage due to insufficient battery power to complete the program.

Note

- Press and hold the Power/Stop button [

Operation

Chair control via OSIM Well-Being App

Full control of the functions is available via OSIM Well-Being App.

1. Download the OSIM Well-Being App, and follow the registration steps.
2. After registration, press [], to exit Sleep mode.
3. Locate the product's Bluetooth ID within the App for pairing, the Bluetooth LED indicator will emit a Solid Blue Light once connected.
4. Follow the in-App instructions to pair with the product.
5. Once completed, follow the in-App instructions to set-up your OSIM account and use the App functions.

Note

- Please keep OSIM Well-Being App running at the foreground at all times when massaging.

App download


OSIM Well-Being App download

Download the OSIM Well-Being App to enjoy additional product features exclusively through the App.

Method 1 :

a: You can directly open OSIM Well-Being App via NFC tag.

Tap your smart device against the NFC tag on the product and it will open up the OSIM Well-Being App.

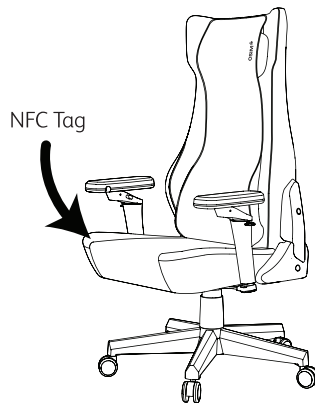b: You will be redirected to the OSIM Well-Being App if it is installed on your smart device.

c: If the OSIM Well-Being App is not installed on your smart device, you will be redirected to the App Store, Google Play or the OSIM Well-Being App introduction page (depending on your mobile operating system) for easy download.

Method 2:


Scan the QR code below to download and install App for free.

App download

Method 3:

Go to the respective App stores (depending on your operating system and region), search for “OSIM Well-Being” to download and install App.

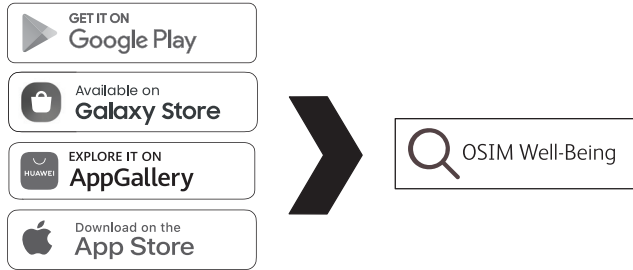Note

- Please check the NFC location on your smart device before using this feature as it varies for different models.
- You may need to remove the cover on your smart device for better sensitivity.

Functionality of OSIM Well-Being App

- You are able to control the product through the OSIM Well-Being App.
- New features will be highlighted and made available in the App, so that the product can stay up-to-date with the latest features.

Compatibility of OSIM Well-Being App

- Due to variations in operating systems, non-compatibility could occur. We recommend updating your mobile operating system whenever a new version is available to enjoy the new App features. The minimum requirements for mobile operating system are subjected to change without notice.
- All smart devices require a compatible wireless connectivity (Wi-Fi or cellular data networks).
- For updates on OSIM Well-Being App, please visit the respective App stores.
- Notifications of App updates are encouraged to be enabled, and it is recommended that you install all updates to ensure that the App will function correctly and is kept up to date.

Product care and storage

Chair placement and storage

Environmental factors can strongly affect the product physically and mechanically. Do not expose the product to direct sunlight or in places with extreme heat or air-conditioning vents. These can cause the upholstery to dry up, crack or the color to fade.

If the product is not to be used for a long period of time, cover it with a dust cover and disconnect the power plug from the electrical outlet.

Ensure enough space between the chair and surrounding furniture items and wall for unobstructed movements during chair usage. Direct contact between the exterior of the chair with other furniture items and walls will cause abrasive damage to the upholstery.

Clean the product using a dry lint-free cloth

Before cleaning, ensure that the main switch is turned off and the power plug is disconnected from the electrical outlet.

The upholstery of the product is customarily exposed to the natural body oils of users. The main unit structure and fabric can be cleaned thoroughly with a piece of soft, dry cloth. The common micro-fiber cleaning materials or gloves are very well suited for this routine maintenance process. Such materials and gloves are available for sale in most DIY shops.

Do NOT use wet elements for cleaning

The product is made up of electrical components and mechanical parts, which are all very sensitive to moisture. Under no circumstances should they be exposed to any form of water or moisture. If you use water (eg. a wet cloth) to clean the upholstery, the coating may be eroded and the protection for the product will be diminished.

As for the control panel, it may result in malfunction if you use a wet cloth to clean.

Do NOT use cleaning solvents

Under no circumstances should you use any cleaning agents, like benzene or thinner to clean the product. This may cause the color of the upholstery to fade.

Do not wash the upholstery using washing machine. This may result in the deterioration of upholstery quality.

Do NOT attempt to repair your product

Under no circumstances should you attempt to repair the product, for safety reason. If you suspect a malfunction, always call OSIM Service Center* for further assistance and information.

Troubleshooting

If you encounter problems with the operation of the product, please check the following points to see if they can be rectified. If the product still malfunctions or does not work, please switch off the product and remove the power source connector from the battery pack. Contact OSIM Service Center* for servicing.

Do not disassemble or attempt to repair the product yourself. Doing so will void the warranty and OSIM will not be responsible for any subsequent consequences.

Problem	Possible causes/actions
• Motor sound during operation.	• These sounds originate from the internal mechanism within the product during operation. • This may occur during normal usage and is not due to malfunction.
• The control panel does not function correctly.	• Check the power source connector is secured to the battery. • Ensure that control panel is lighted up. • Ensure that the product is not connected to the OSIM Well-Being App.
• The battery gets abnormally hot.	• Stop using the product immediately, unplug the battery. • Remove battery from chair and place it somewhere without any combustible material. • Contact OSIM Service Center* for assistance immediately.
• The massage stop midway thru a 15 minutes program.	• The auto timer of 15 minutes has lapsed. Press any program button to restart the massage if the control panel is lit up. • If the product fails to restart, and control panel is off, unplug battery for charging (or replace with your fully charged spare battery). If control panel lit up, select auto program again. • If the same problem persists, unplug the battery and contact OSIM Service Centre*.

Troubleshooting

Problem	Possible causes/actions
• The product is powered up, but cannot operate, the control panel LED are blinking in certain sequence / pattern.	• Memorize or do a recording of the sequence of the LED patterns. Unplug the battery and remove from chair. • Contact OSIM Service Centre* for assistance.
• If you encounter any problems other than the ones stated above, or the same problem persists.	• Unplug the battery and remove from chair and contact OSIM Service Centre* for assistance.

* Please visit us at www.OSIM.com for information on the nearest service center.
Product warranty is valid only in the country of purchase.

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© 2025 All rights reserved. www.OSIM.com
Designed & Developed in Singapore by OSIM International Pte. Ltd.